











ISO 14001 REGISTERED

QMS

4Front | Features and Options



Smooth Ride Suspension

Superior, automotive-grade 'SRS' suspension (Smooth Ride Suspension) CASE (Caster Angle Sensor Encoder) technology

70Ah Batteries

More miles per charge! 70 Ah batteries allow for increased range, for the most active user. Easy access through front shroud.

Front-Wheel Drive

Quiet and responsive the front-wheel drive system of the 4Front features exceptional drive performance.

Turning Radius

Exceptional tight space turning radius of 24.5" (62.23 cm).

Motor Package

2-Pole advanced responsive motors 6 Mph motor and featuring the new Q-Logic 3 control system.

Measurements

14" Drive wheel, 9" low impact solid castor rear wheels Length: 40" (Without legrests) | Width: 24" | Ground Clarence: 3".

Seating Solutions



TRU-Balance® 3 Seating

- 12" x 12" to 22" x 22" seat size options
- Make fine adjustments while client is in the chair
- Accepts a wide variety of after market seats and backs
- Endomorph and mesomorph seating, allows for independent adjustment of seat width and backrest width
- Premium accessories available

Controls



Q-Logic 3 Controls

- Dongle Free Programming
- Bluetooth included
- SD card slot for programming, firmware updates, and storing images.
- New simple easy to use enhanced display
- Smart easy to use interface

4Front | Features and Options

Power Options



TRU-Balance® 3 Power Tilt

- 0-50° of tilt (Anterior tilt available)
- · Lower standard seat-to-floor height
- Easy adjustment and durability
- Maintains centre of gravity shift

TRU-Balance® 3 Power Recline

- 85° 175° of recline
- · Adjustable shear reduction
- Compatible with after-market backrest interface

Powered Elevating leg rests

- Articulating foot platform
- · Articulating swing-away leg rests
- Can be operated in conjunction with TB3 recline



10" Safe Seat Elevation

- Safe seat elevation allows 10" of seat elevation while at a speed of 3.2 mph.
- Independent wheel and castor travel, provides additional stabilisation.

Additional Features



TRU-Balance 3 Cantilever Flip-up armrests

- Quick height adjustment
- Can be internally or externally rotated for optimal positioning
- Multiple armrest pad options



TRU-Balance® 3 Centre mount foot platform

- · Fixed or powered articulating
- Modular footplate extension kits
- Multi-adjustable calf padsPaediatric & adult footplate options



Swing-away legrests

- 60° 90° standard swing-away
 Paediatric & adult footplate
- options
 Powered Articulating
- · Manual elevating



Attendant/Dual control options

- Full proportional joystick
- Quick and easy to use control interface
- Control the chair at the attendants speed



Body positioning

- A full range of multi adjustable swing-away laterals, offering advanced positioning of the trunk
- Designed to easily allow head to toe positioning using any of the Stealth positioning accessories



After-market backrest interface

- Allows for an unprecedented level of compatibility between TRU-Balance 3 back rest and most after-market backrests
- Retains full functionality, including shear reduction



Premium accessories

- Transfer bars
- · Mobile phone holder
- USB Port (available with Q-Logic 3)
- Clothes guard
- · Additional accessories available



Full range of control devices available

- Switches
- Sensors
- Proportionals
- i-DriveTM



4Front | Colour Options



Pretty in Pink



Red Corvette



Lime Light



Midnight Blue



Yellow Blaze



White Satin





Night Train





Silver







Grey Street




Military Green

Technical Specifications







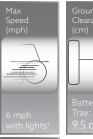



















Warranty

- Lifetime limited warranty on frame.
- 1-year limited warranty on electronics.
- 18-month warranty on drive control system.
- 2-year limited warranty on drivetrain.

Notes

- Speed varies with user weight, terrain type, battery charge, battery condition and tyre pressure.
- Due to manufacturing tolerances and continual product improvement, this specification can be subject to variance of + or - 3 %. Dimensions listed are for power base only. Overall measurements will vary based on seating and accessory selections.
- Battery weight may vary based on manufacturer.

Get Social with Quantum



FOLLOW THE LIVES OF OUR QUANTUM STARS, AND GET THEIR PERSPECTIVE OF THE Q6 EDGE 2.0

www.quantumrehab.co.uk/quantum stars



Ensure you have the most up to date information about Quantum Rehab by subscribing to our social media pages

#ileveluk

#myquantumproduct

#socialmobility

Further information can obtained through our web site, or prescription forms. All information is subject to change without notification. Please consult Quantum with any queries you may have.